

0366

REPORT DOCUMENTATION PAGE

Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing this collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden to Department of Defense, Washington Headquarters Services, Directorate for Information Operations and Reports (0704-0188), 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to any penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number. PLEASE DO NOT RETURN YOUR FORM TO THE ABOVE ADDRESS.

1. REPORT DATE (DD-MM-YYYY) 05 May 00		2. REPORT TYPE Final Technical Report		3. DATES COVERED (From - To) 05 Aug 96 - 04 Feb 00	
4. TITLE AND SUBTITLE Infrared Cavity Ringdown Laser Absorption Spectroscopy: Metal-containing Clusters and HEDM Molecules				5a. CONTRACT NUMBER	
				5b. GRANT NUMBER =F49620-96-1-0411	
				5c. PROGRAM ELEMENT NUMBER	
6. AUTHOR(S) Saykally, Richard J.				5d. PROJECT NUMBER	
				5e. TASK NUMBER	
				5f. WORK UNIT NUMBER	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) The Regents of the University of California University of California C/o Sponsored Projects Office 336 Sproul Hall #5940 Berkeley, CA 94720-5940				8. PERFORMING ORGANIZATION REPORT NUMBER	
9. SPONSORING / MONITORING AGENCY NAME(S) AND ADDRESS(ES) AFOSR/NL 810 N Randolph St Room 732 Arlington, VA 22203-1977				10. SPONSOR/MONITOR'S ACRONYM(S) AFOSR/NL	
				11. SPONSOR/MONITOR'S REPORT NUMBER(S)	
12. DISTRIBUTION / AVAILABILITY STATEMENT Approved for public release; distribution is unlimited.					
13. SUPPLEMENTARY NOTES					
14. ABSTRACT Efforts focused on four activities related to the central objective of studying metal-containing clusters, HEDM molecules, and clusters of ceramic materials. Nineteen AFOSR supported papers were published on the following topics (including five journal cover articles): 1) Development of high resolution (100 MHz), high sensitivity pulsed infrared cavity ringdown spectroscopy as a general tool for trace species characterization. 2) Infrared cavity ringdown system with medium (900 MHz) resolution constructed and operational. 3) IR studies of hydrogen bonding in water, methanol, ethanol, and butanol clusters carried out. 4) More difficult IR experiments on laser vaporized carbon and silicon carbide initiated.					
15. SUBJECT TERMS Laser Spectroscopy, Ceramic Materials, Clusters.					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT UU	18. NUMBER OF PAGES 13	19a. NAME OF RESPONSIBLE PERSON Richard J. Saykally
a. REPORT U	b. ABSTRACT U	c. THIS PAGE U			19b. TELEPHONE NUMBER (include area code) (510) 642-8269

Standard Form 298 (Rev. 8-98)
Prescribed by ANSI Std. Z39.18

FINAL TECHNICAL REPORT

AFOSR GRANT #F49620-96-1-0411

**"Infrared Cavity Ringdown Laser Absorption Spectroscopy: Metal-containing
Clusters and HEDM Molecules"**

(Period Covered: 08/05/96 – 02/04/00)

Richard J. Saykally - Principal Investigator

Department of Chemistry
University of California, Berkeley
Berkeley, CA 94720-1460

Efforts focused on four activities related to the central objective of studying metal-containing clusters, HEDM molecules, and clusters of ceramic materials. Nineteen AFOSR-supported papers were published on the following topics (including five journal cover articles): 1) Development of high resolution (100 MHz), high sensitivity pulsed infrared cavity ringdown spectroscopy as a general tool for trace species characterization. 2) Infrared cavity ringdown system with medium (900 MHz) resolution constructed and operational. 3) IR studies of hydrogen bonding in water, methanol, ethanol, and butanol clusters carried out. 4) More difficult IR experiments on laser vaporized carbon and silicon carbide initiated.

Ph.D. Degrees awarded to:

Dr. Joshua B. Paul (Fall, 1998)

Dr. Robert A. Provencal (Fall, 1999)

M.S. Degree awarded to Chris Chapo (Spring, 1999)

Technical Summary

Infrared spectroscopy is a powerful and time-honored technique for characterizing chemical materials. State-of-the-art technology for IR spectroscopy employs Fourier transform spectrometers, color center, diode, and difference frequency lasers, and pulsed optical parametric amplifier (OPO/OPA) systems. These widely used approaches all suffer from either complexity, restricted tuning range, or relatively low sensitivity. IR cavity ringdown laser absorption spectroscopy (CRLAS) with pulsed laser sources promises an impressive combination of convenience, high sensitivity, and very broad tuning range. The major focus of our work during this grant period has been to advance the technique of IR-CRLAS, and to apply it to a variety of trace species produced in pulsed supersonic beams, including metal-containing molecules, hydrogen bonded clusters, and covalent carbon and silicon carbide clusters.

A medium (ca. 900 MHz) resolution pulsed IR CRLAS spectrometer employing Stimulated Raman Scattering (SRS) of a pulsed dye laser with a continuous tuning range from visible to 10 μm , was constructed and used in detailed studies of water methanol, ethanol, and butanol clusters. These studies were extended to the more difficult class of laser-vaporized carbon and silicon carbide clusters, with the goal of applying these to ceramic materials and HEDM species. The HEDM molecule Al NH was also studied.

A high resolution (100 MHz) version of this experiment was designed built, and tested. It features an external cavity diode laser seeded Alexandrite laser to produce 40 MHz pulses in the 700-800 nm region and SRS to downconvert into the IR, and has tunability from far IR to the visible. Numerous invited cover articles were published on these systems and their initial applications. Initial studies of metal-containing molecules (WO and Cr_2) were carried out collaboratively with the Bondybey group at the Technical University in Munich in the visible region of the spectrum.

A Medium-Resolution (900 MHz) IR CRLAS Spectrometer

The design of this system is shown in Figure 1, and has been described in numerous feature articles [155, 181, 190, 191, 204, 220]. It employs a commercial YAG-pumped dye laser at visible wavelengths and a long multi-pass high pressure hydrogen cell for stimulated Stokes scattering downconversion into the IR. This system produces ca. 5 ns pulses with 100-500 μJ of energy over the range 10 μm – 2 μm . With available "supermirrors," we achieve fractional absorption sensitivity of $\sim 1 \times 10^{-6}$ throughout most of this region.

We have first exploited this system for the study of water and alcohol clusters in pulsed planar supersonic expansion. After perfecting the technology, we were able to make the first direct absorption measurements of H_2O and D_2O dimers, trimer, tetramer, pentamer, and hexamers.

This allowed us to exact absolute number densities for the clusters in the expansion and to assess the capability of this method as an *in situ* probe of water clusters and other species in equilibrium environments. We were also able to carry out the first systematic study of alcohol (methanol, ethanol, and butanol) clusters in the mid-IR as a result of combining the H₂ SRS source with our CRLAS detection scheme (Figure 2).

The next step was to perform the much more difficult IR-CRLAS studies of clusters generated by laser vaporization in supersonic beams. In order to develop this capability, we carried out studies of the pure carbon clusters C₃, C₅, and C₉ previously studied by our diode laser experiments in order to perfect the combination of these laser plasma sources with ringdown detection. We seek to exploit the ultra-high sensitivity and rapid scan capability of CRLAS, which will be indispensable for the proposed work with ceramic clusters. The ν_3 band of C₅ is shown in Figure 3 for comparison.

We are presently nearing the ability to begin study of metal carbide clusters and HEDM species. As an initial step we initiated a study of Al NH, itself of possible HEDM interest and previously observed by IR spectroscopy in a cryogenic matrix. Representative spectra, measured in a plasma of laser vaporized aluminum with NH₃ in the He buffer gas, are shown in Figure 4. Analysis is still in progress.

Over the last 1 1/2 years, we have devoted much effort to the design, construction, and testing of a high resolution IR-CRLAS experiment. We decided to base the design on a tunable Alexandrite laser (Light Age Inc.) which provides >150 mJ/pulse with a 50 MHz transform limited bandwidth. This high spectral brightness is converted into the IR by SRS in a combination of H₂ and D₂ Raman cells. This scheme will allow us to cover the entire mid-IR with narrowband pulsed tunable radiation, permitting us to measure infrared spectra of virtually all molecules with Doppler-limited resolution (a 10x improvement over the present system). This new CRLAS system is currently operational with narrow linewidths over a limited wavelength range. We are working to extend the range of this high resolution system.

A High Resolution (100 MHz) IR CRLAS Spectrometer

Most of our efforts over the last two years have been devoted to the design, construction, and testing of a high resolution IR-CRLAS experiment. This would achieve a 10-30 fold increase in sensitivity and provide Doppler-limited spectra, like diode lasers. We decided to base the design on a tunable injection-seeded Alexandrite laser (Light Age Inc.), rather than on the OPO/OPA system originally requested in this proposal. We chose this route because the combination of the Alexandrite laser and SRS I H₂ and D₂ cells promised continuous spectral coverage from far IR (ca. 500 cm⁻¹) to the UV, with narrow linewidths – probably limited by the SRS downconversion process. Both the tuning range and the spectral resolution can thus greatly surpass those available

with commercial OPO/OPA systems. However, this approach required great care in the optical design and the detailed exploration of the SRS process with regard to linewidth limitations. The design of the system is shown in Figure 5.

We have now demonstrated (Figure 6) that this system does indeed produce tunable IR pulses with linewidths less than 250 MHz. Moreover, it appears that ca. 500 μ J pulses can be produced over much of the IR (measured for 1900-2000 cm^{-1} but currently limited by the availability of the diode seed lasers). We are just now beginning experiments on laser vaporized species with this system to test the synchronization of the laser pulses with the pulsed production of the sample, which has also proven to be a difficult modification of the Alexandrite system, achieved by Light Age Inc. We have worked closely with them in the design of this spectrometer, which should be producing spectra of new species very soon.

Molecules Studied

A. Metal-Containing Species

A major object of the technical development described above is the study of metal-containing species, particularly clusters comprising ceramic materials. We have collaborated with the Bondybey group in Munich in the study of simple high temperature metal oxides (e.g. tungsten oxide) [218], by CRLAS in the visible, as they have developed sophisticated technology for producing such refractory species in pulsed supersonic jets. Chromium dimers were also studied, but attempts to extend this approach to titanium have not yet worked [222]. The plan is to adapt these sources for use in our IR cavity ringdown spectrometers.

B. HEDM Species

The molecule Al NH has been proposed as a possible HEDM molecule, being observed previously in a cryogenic matrix. We have measured rotationally resolved spectra near 1890 cm^{-1} , as shown in Figure 4, by laser vaporizing aluminum with NH_3 in the helium buffer gas. Carbon clusters and silicon carbide clusters have also been studied by IR CRLAS at "medium resolution" (Figure 3).

C. Hydrogen Bonded Clusters

As discussed above, a series of studies of water clusters and alcohol clusters, were carried out with the goal of characterizing the differences in hydrogen bonding that are effected as the H in a molecule is replaced by sequentially larger hydrophobic groups. We decided to carry out this study largely to gain experience in the difficult technology of the IR CRLAS system before undertaking the much more demanding laser vaporization experiments.

Publications Supported by This Grant
(numbers correspond to complete publication list)

- 155.* J.J. Scherer, J.B. Paul, A. O'Keefe, and R.J. Saykally, "Cavity Ringdown Laser Absorption Spectroscopy: History, Development, and Application to Pulsed Molecular Beams," *Chemical Reviews* **97**, 1, 25-51 (1997).

***Cover Article.**

189. J.B. Paul, C.P. Collier, J.J. Scherer, A. O'Keefe, and R.J. Saykally, "Direct Measurement of Water Cluster Concentrations by Infrared Cavity Ringdown Laser Absorption Spectroscopy," *J. Phys. Chem. A* **101**, 5211 (1997).

- 190.* J.B. Paul and R.J. Saykally, "Cavity Ringdown Laser Absorption Spectroscopy," *Analytical Chemistry A* **69**, 287 (1997).

***Cover Article.**

- 191.* J.B. Paul, J.J. Scherer, A. O'Keefe, and R.J. Saykally, "Cavity Ringdown Laser Absorption Spectroscopy," *Laser Focus World* **33**, 71-80 (1997).

***Cover Article.**

196. A.I. Boldyrev, J. Simons, J.J. Scherer, J.B. Paul, C.P. Collier, and R.J. Saykally, "On The Ground Electronic States Of Copper Silicide and Its Ions," *J. Chem. Phys.* **108**, 5728 (1998).

- 204.* J.B. Paul, R.A. Provencal, C. Chapo, E. Michael, A. Pettersson, and R.J. Saykally, "Infrared Cavity Ringdown Laser Absorption Spectroscopy of Transient Species in Pulsed Supersonic Expansions," *1997 ACS Symposium Series 720 on Cavity-ringdown spectroscopy: an ultratrace-absorption measurement technique*, Kenneth W. and Marianna A. Busch, editor, pp. 162-173 (1999).

***Cover Article.**

206. J.B. Paul, R.A. Provencal, and R.J. Saykally, "Characterization of the (D₂O)₂ Hydrogen-Bond-Acceptor Antisymmetric Stretch by IR-Cavity Ringdown Laser Absorption Spectroscopy," *J. Phys. Chem. A* **102**, 3279 (1998).

208. R.A. Provencal, J.B. Paul, E. Michael, and R.J. Saykally, "Cavity Ringdown Laser Absorption Spectroscopy," *Photonics Spectra* **32**, 159-166 (1998).

209. J.B. Paul, R.A. Provencal, C. Chapo, A. Pettersson, and R.J. Saykally, "Infrared Cavity Ringdown Spectroscopy of Water Clusters: O-D Stretching Bands," *J. Chem. Phys.* **109**, 10201 (1998).

213. J.B. Paul, R.A. Provencal, C. Chapo, K. Roth, R. Casaes, and R.J. Saykally, "Infrared Cavity Ringdown Spectroscopy of the Water Cluster Bending Vibrations," *J. Phys. Chem. A* **103**, 2972 (1999).

214. C.J. Chapo, J.B. Paul, R.A. Provencal, K. Roth, and R.J. Saykally, "Is Arginine Zwitterionic or Neutral in the Gas-Phase? Results from IR Cavity Ringdown Spectroscopy," *J. Am. Chem. Soc.* **120**, 12956-12957 (1998).

217. R.A. Provencal, J.B. Paul, K. Roth, C. Chapo, R.N. Casaes, R.J. Saykally, G.S. Tschumper, and H.F. Schaefer, III, "Infrared cavity ringdown spectroscopy of methanol clusters: Single donor hydrogen bonding," *J. Chem. Phys.* **110**, 4258 (1999).

218. D. Kraus, R.J. Saykally, and V.E. Bondybey, "Cavity Ringdown Laser Absorption Spectra of Tungsten Oxide," *Chem. Phys. Lett.* **295**, 285-288 (1998).

220.* R.A. Provencal, J.B. Paul, C.N. Chapo, and R.J. Saykally, "Cavity Ringdown Laser Absorption Spectroscopy," *Spectroscopy* **14**, 24 (1999).

***Cover Article.**

221. A. O'Keefe, J.J. Scherer, J.B. Paul, and R.J. Saykally, "Cavity-Ringdown Laser Spectroscopy: History, Development, and Applications," *1997 ACS Symposium Series 720 on Cavity-ringdown spectroscopy: an ultratrace-absorption measurement technique*, Kenneth W. and Marianna A. Busch, editor, pp. 71-92 (1999).

222. D. Kraus, R.J. Saykally, and V.E. Bondybey, "Cavity ringdown spectroscopy search for transition metal dimers," *Chemical Physics* **247**, 431-434 (1999).

225. R.A. Provencal, K. Roth, J.B. Paul, C.N. Chapo, R.N. Casaes, R.J. Saykally, G.S. Tschumper, and H.F. Schaefer, III, "Hydrogen Bonding in Alcohol Clusters: A Comparative Study by Infrared Cavity Ringdown Laser Absorption Spectroscopy," *J. Phys. Chem. A* **104**, 1423 (2000).

229. D.R. Wagner, and R.J. Saykally, "Ringing the changes," *Chemistry in Britain* **36**, 47-49 (July 2000).

Submitted/In Press Papers

M00.4 E. Michael, C.J. Keoshian, D.R. Wagner, S.K. Anderson, A.J. Huneycutt, and R.J. Saykally, "Infrared Water Lasers," *Nature* (submitted 2000).

Figure Captions

Figure 1. Design of Medium Resolution IR Cavity Ringdown Spectrometer

Tunable infrared radiation is generated by Raman shifting an ND:YAG pumped dye laser. The tunable infrared light is injected into the ringdown optical cavity where the beam intersects a supersonic expansion in which the transient molecules are formed. Light exiting the cavity is monitored using a detector (InSb or HgCdTe). The signal is amplified, digitized, and sent to a PC computer for determination of cavity losses.

Figure 2. CRLAS Spectra of Methanol Clusters in the O-H Stretching Region

Spectrum obtained using the medium resolution Cavity Ringdown spectrometer. Features from the dimer to tetramer were assigned along with a broad band from the pentamer and larger clusters.

Figure 3. The ν_3 Band of C_5

The R(12) to R(26) transitions of the ν_3 band. The sensitivity of the medium resolution CRLAS spectrometer is as good as the previous diode laser experiment. The high resolution CRLAS spectrometer should be 10 times more sensitive.

Figure 4. IR Spectrum of Al NH Molecule

Spectra assigned to Al NH on the basis of IR matrix data. This HEDM molecule was generated by laser vaporization of Al into a mixture of He and NH_3 .

Figure 5. The High Resolution IR Cavity Ringdown Spectrometer

500 micro joules of infrared radiation (20700 cm^{-1} to 1900 cm^{-1}) are generated by stimulated Raman shifting of the Alexandrite laser (100 mJ/40 MHz) in deuterium and hydrogen.

Figure 6. Spectrum of OCS Confirming Instrument Linewidth <250 MHz

Linewidth vary from 350 MHz to 150 MHz. This is probably due to variations in the Alexandrite linewidth during the scan. This problem will be fixed once we install a new OEM.

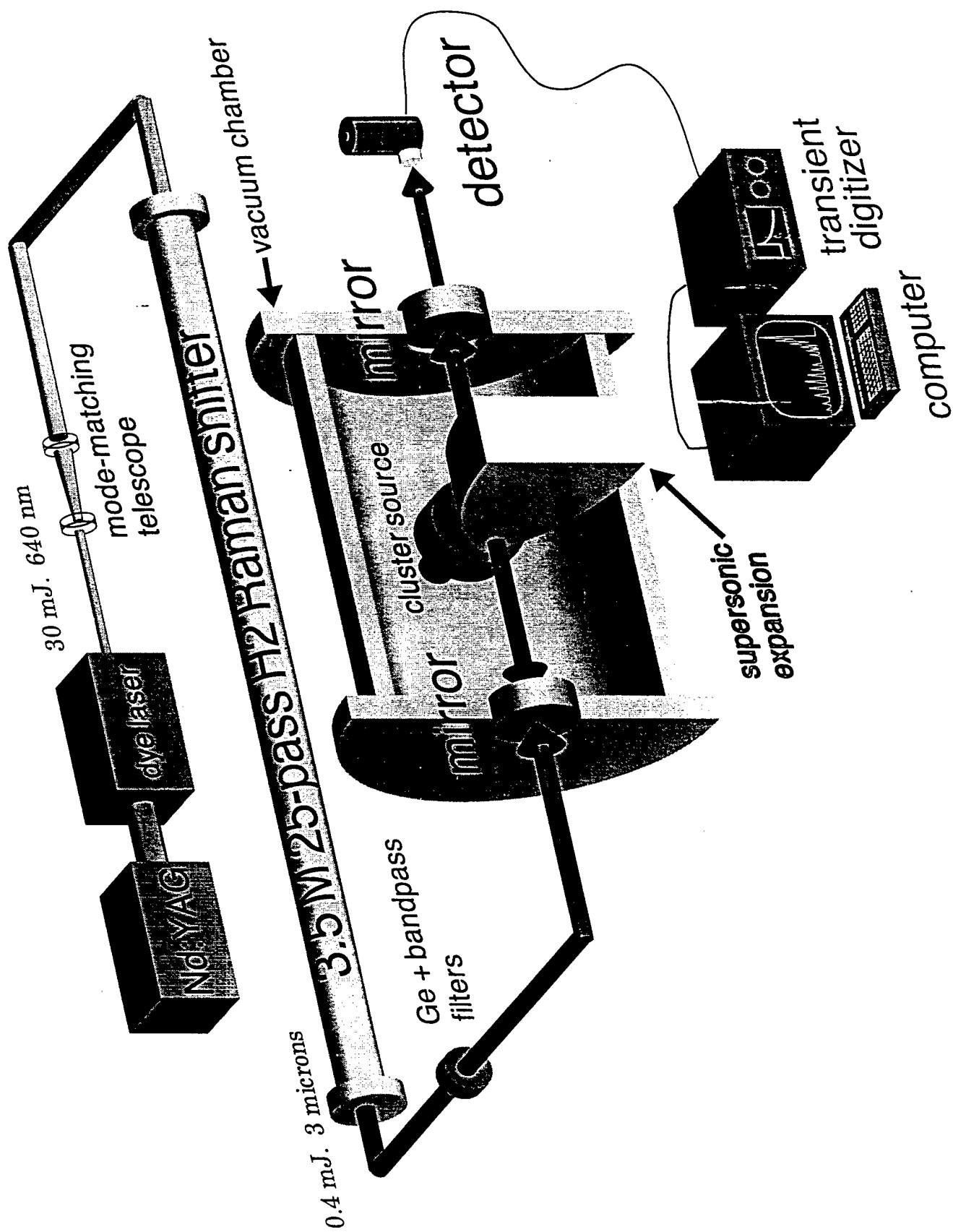


Figure 1

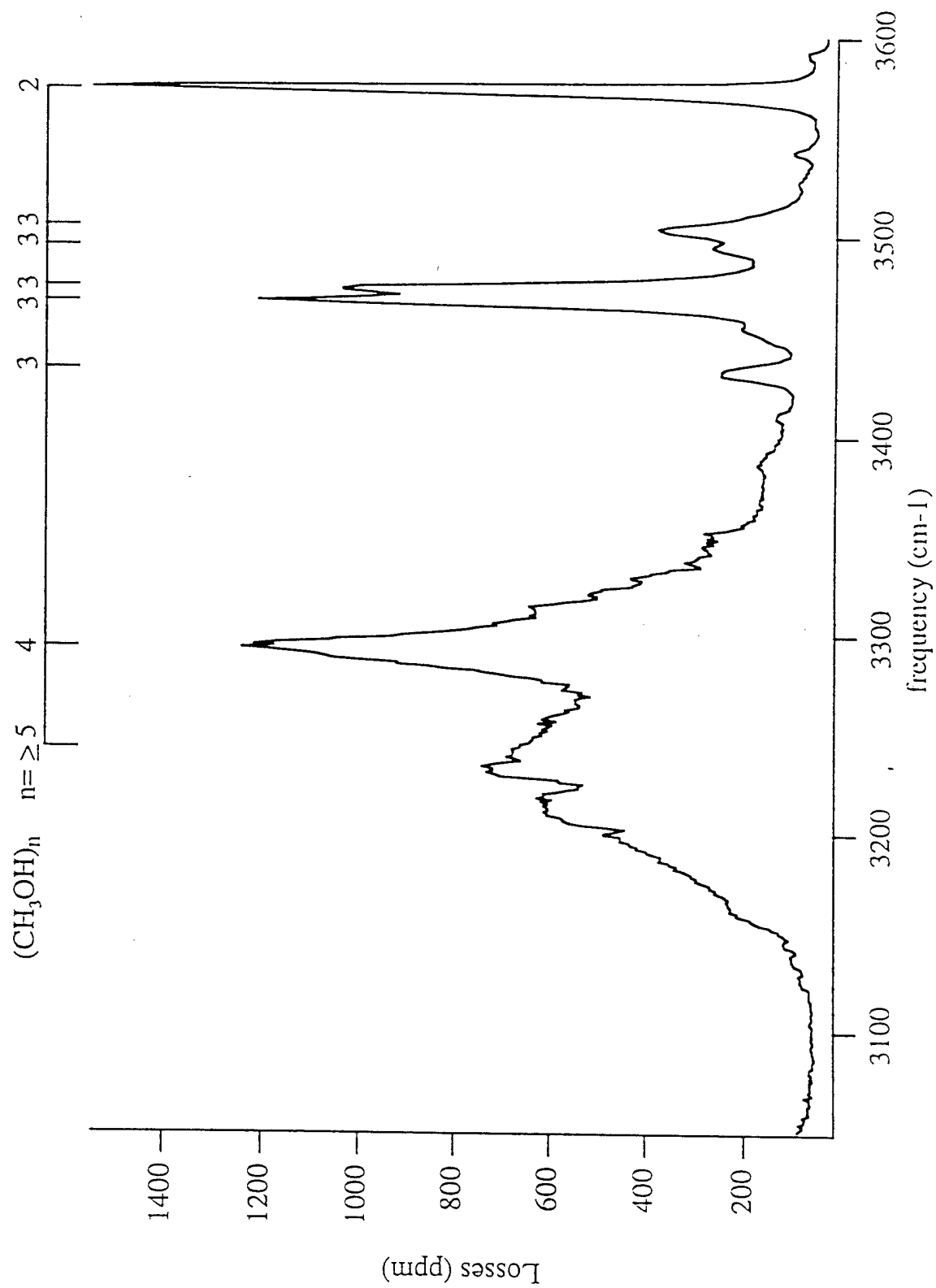


Figure 2

The ν_3 band of C_5

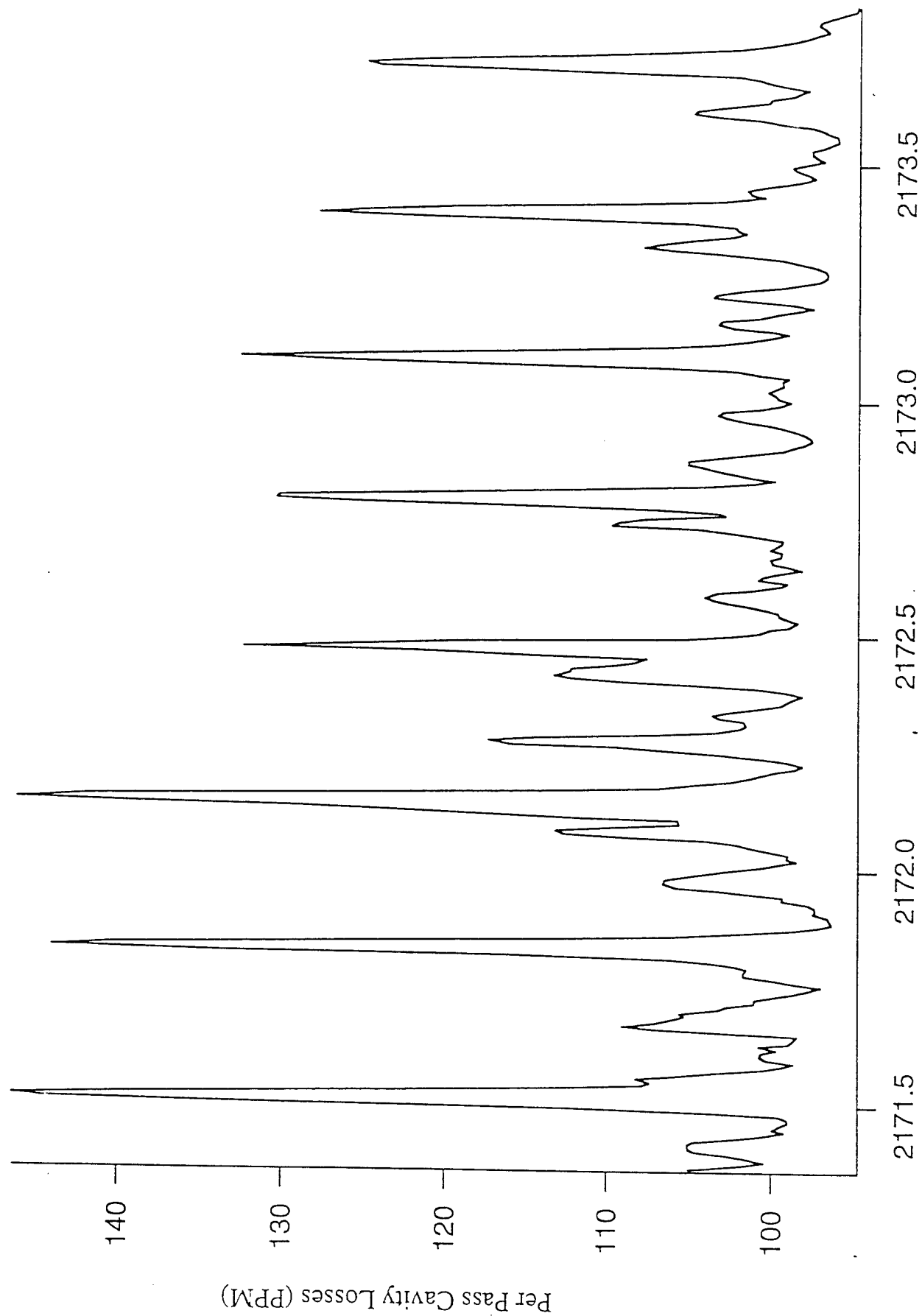


Figure 3

Al-N-H data

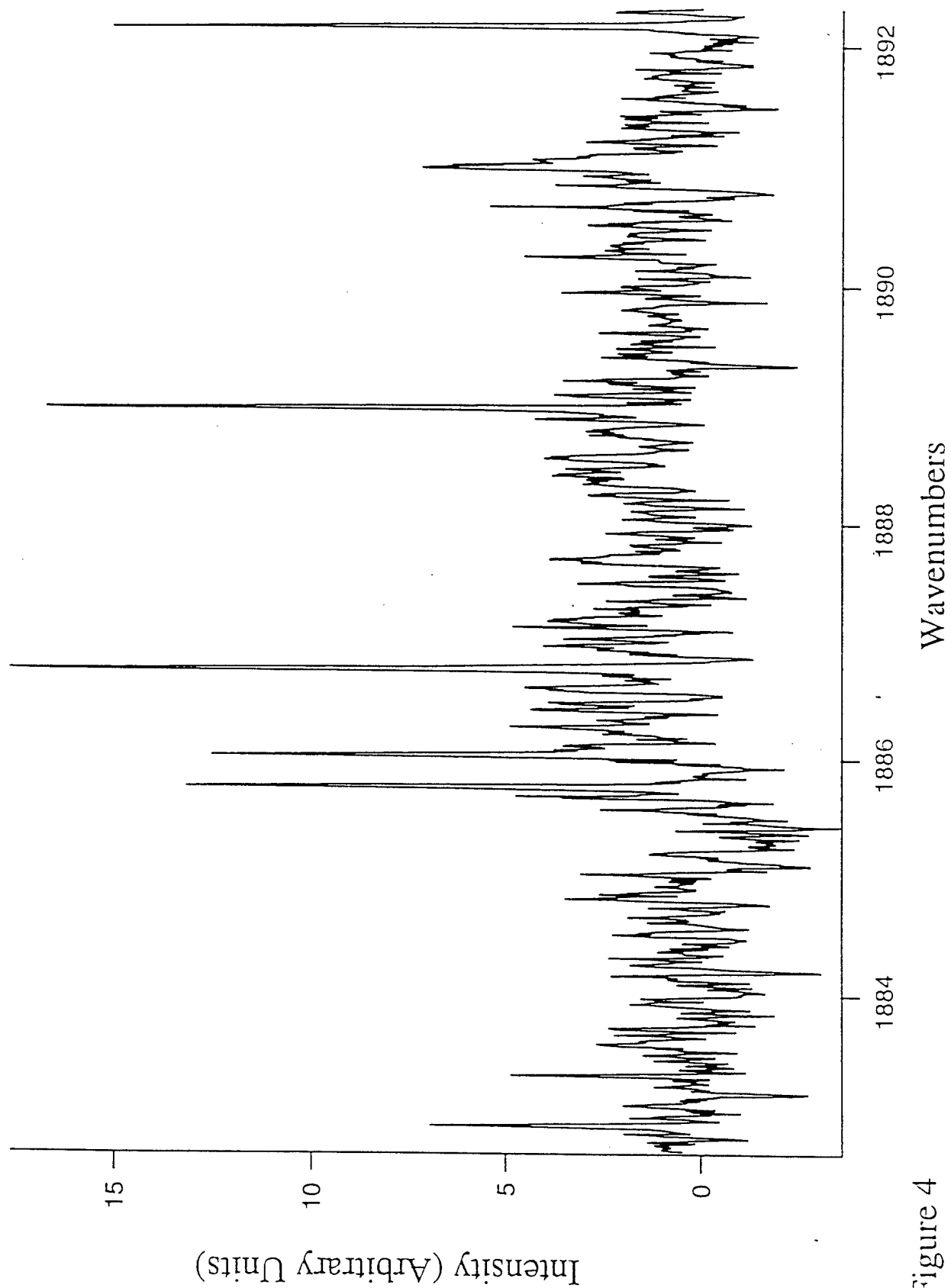


Figure 4

High Resolution Cavity Ringdown Spectrometer

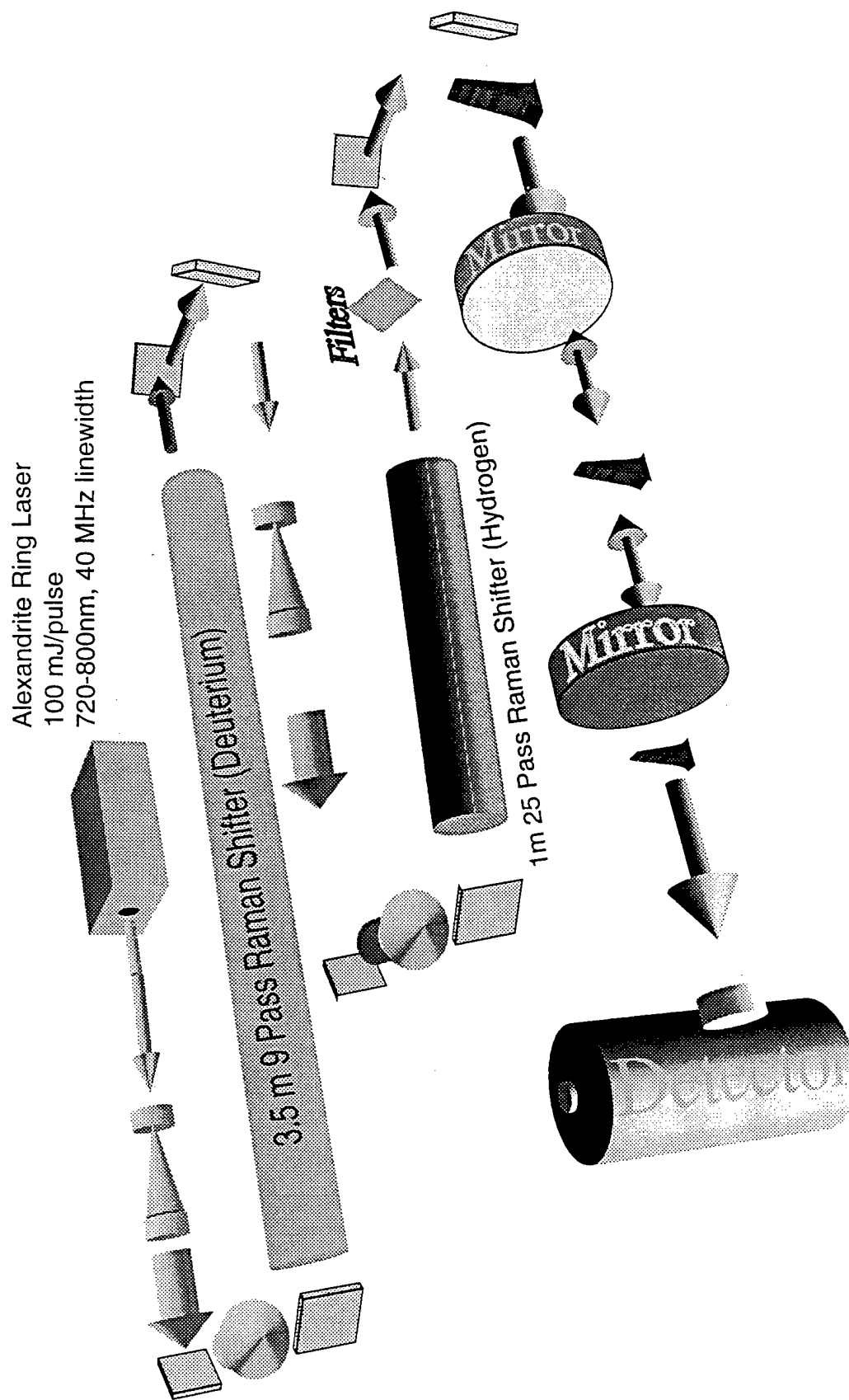


Figure 5

OCS in a Static Cell

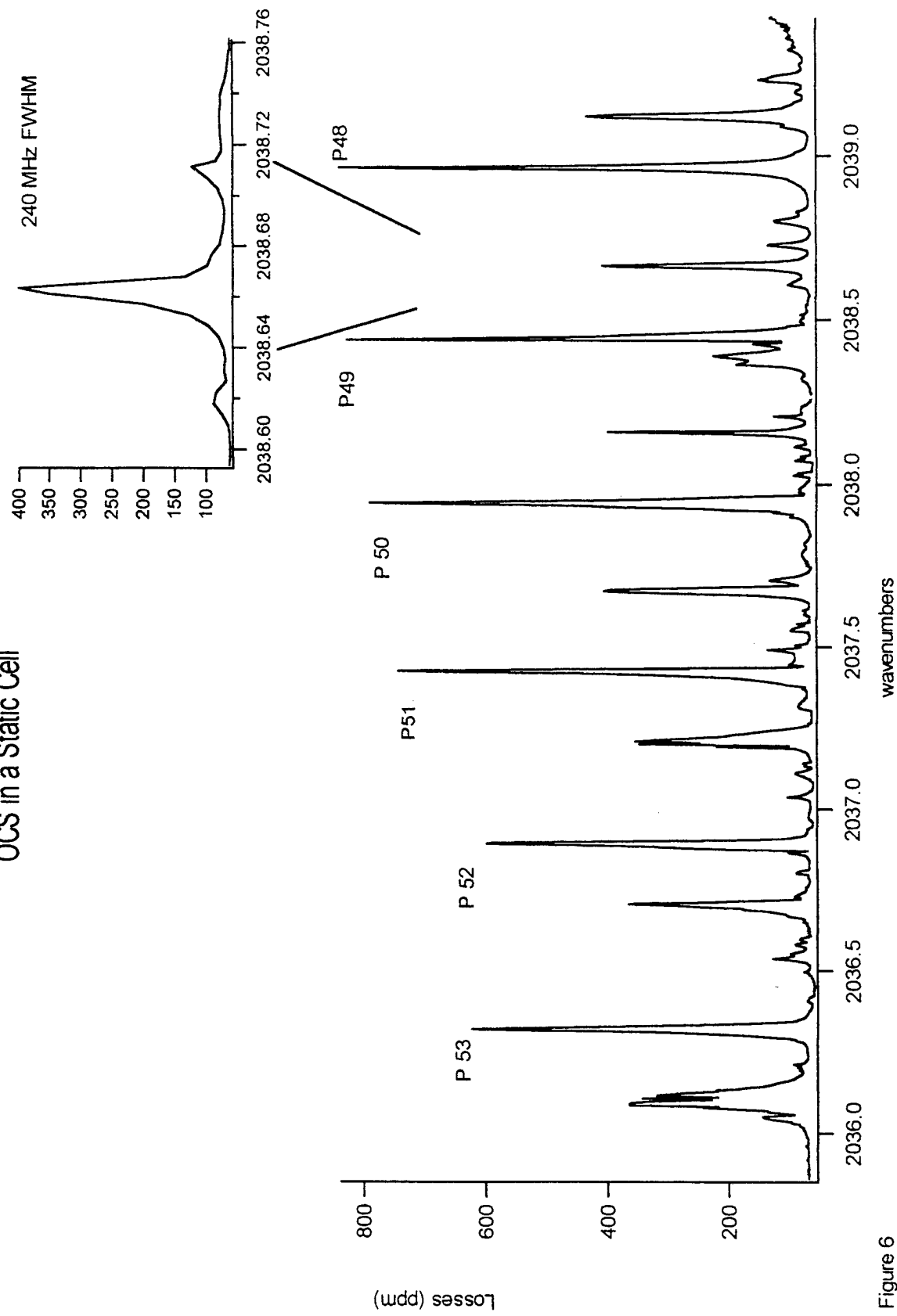


Figure 6